

The Invacare® Top End® Eliminator™ OSR racing chair is designed for speed. The super-stiff roll-formed horizontal maintube, combined with a relaxed headtube angle for exceptional stability and aerodynamics, make it the choice of champion racers like Ernst van Dyk, Krige Schabort, Kurt Furnley, Chantal Petitclerc and Edith Hunkeler. So if you have your sights set on winning, the Eliminator OSR racing chair will help you get to the finish line *first*.



Options

- Carbon fiber rear wheels
- Corima™ carbon fiber front wheel
- Alignment gauge
- Wheel bag
- · Racing suit
- Solid seat/kneel plate
- Fixed footrest
- Cushion
- Cordless speedometer
- Helmet
- 18" front wheel



Oversized rake for maximum speed and cornering.



The precision track compensator system ensures maximum speed on the road and around each turn of the track.



Features

- 6061 T6 oval tube aluminum frame
- Integral wrap-around fenders
- 4 cage designs: V, Open-V, U or I-cage
- Kneeling, amputee or traditional upholstery
- 700C or 26" high performance rear wheels
- Caliper hand brake
- Aluminum side panels
- Straight axle tube with rear wheel alignment _____
- Aluminum fork
- Precision track compensator system
- 20" front wheel with tubular tire
- Tire-coated 3/8" diameter aluminum rims
- Aero bars
- Powder coat paint




Chantal Petitclerc, Team Invacare. The most medaled track and field athlete in history.





Team Invacare's Edith Hunkeler streaks to set a new course record at the New York City Marathon®.



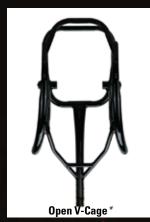

* Shown with cut out camber tube option for more legroom





Standard vs. Kneel

The Invacare® Eliminator™OSR racing chair is custom built to make the most of each athlete's potential by matching up the athlete's ability with the optimum position. The standard position is recommended for athletes who have limited balance and need more of a bucket or squeeze position. The athlete can sit in the traditional position with his/her feet on a footplate or use it in the kneeling position which makes it the ideal choice for beginners. The kneeling position is best for experienced athletes with good balance who only kneel in their chair or need an amputee seat.

Seat Cage Styles

Invacare Top End offers different cage styles to maximize each athlete's transfer skills while maintaining frame stiffness. The I-,V-and Open V-cage styles are offered on the standard OSR racing chair. The Open V-cage style is the most popular because it allows more leg room and easier transfers. The U-cage requires a Top End specialist to fit it properly to the racer. The U-, V- and Open V-cage styles are available on the kneel-only OSR racing chair.

Top End

Sports and Recreation Products 4501 63rd Circle North Pinellas Park, Florida 33781-5914 (727) 522-8677 (800) 532-8677 Fax: (727) 522-1007

Invacare Corporation www.invacare.com

USA

One Invacare Way Elyria, Ohio 44035-4190 (800) 333-6900

Canada

570 Matheson Blvd. E. Unit 8 Mississauga, Ontario L4Z 4G4 Canada (800) 668-5324

©2010

Invacare Corporation. All rights reserved. Trademarks are identified by the symbols [™] and [®]. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. New York City Marathon is a registered trademark of New York Road Runners Club, Inc. Corima is a registered trademark of Corima, S.A. Zipp is a registered trademark of Compositech, Inc. Specifications are subject to change without notification.

Form No. 01-192 Rev 08/10



